
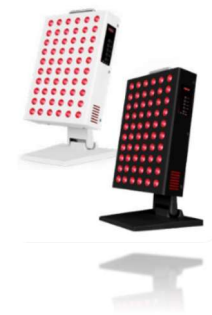


All About Korean Beauty

Sunplus EW300 Plus User Manual

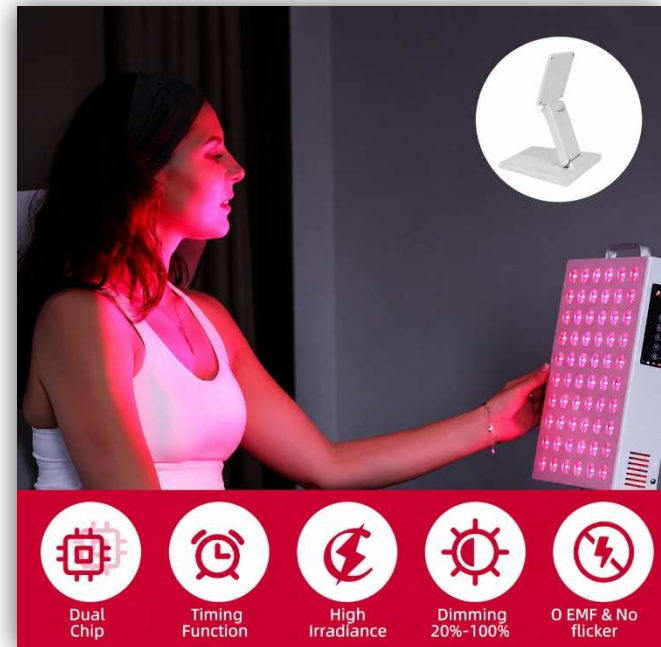

Contents

Welcome..... 1

What’s included 1

Warnings..... 2

Cautions 3

Quick Start Guide 4

How Red-Light Therapy Works..... 5

Key Benefits 5

Technical Specifications 5

Recommended Usage Guidelines..... 6

Maintenance 6

Disclaimer..... 6

Warranty 7

Frequently Asked Questions..... 8

overall wellness benefits. Red light primarily targets the skin’s surface, while near infrared light penetrates deeper into the tissues and muscles.

The main exception is for individuals with melasma or certain pigmentation concerns, as near infrared light may potentially worsen discoloration in some cases. If you have melasma, it is generally recommended to use Red Light Only unless otherwise advised by a healthcare professional.

Q: Can red light therapy be used alongside other skincare products or treatments?

A: Yes, red light therapy is commonly used alongside skincare products and other wellness or aesthetic treatments. However, some ingredients, procedures, or medications may increase sensitivity to light. For best results and safety, we recommend consulting with a healthcare professional or dermatologist to confirm compatibility with your current skincare routine or treatments.

Frequently Asked Questions

Q: How long should I use the SunPlus EW300 Plus for?

A: Treatment time will vary depending on the area being treated and the distance from the device. For detailed dosing recommendations, please refer to the section titled *Recommended Usage Guidelines*. In most cases, sessions of approximately 10 minutes or less are sufficient.

Q: Should I wear eye protection?

A: Eye protection is generally not required when using the SunPlus EW300 Plus. However, it is recommended to avoid staring directly into the LEDs, as the light may feel uncomfortable to some users. Most individuals quickly adjust to the brightness and may comfortably keep their eyes open during treatment, provided no discomfort or sensitivity is experienced.

Q: Why can't I see the infrared LEDs working?

A: The near-infrared (NIR) LEDs on the SunPlus EW300 Plus may appear as though they are not functioning because near-infrared light is outside the visible spectrum and cannot be seen by the human eye. In low-light or dark conditions, you may notice a faint glow from the LEDs, which indicates they are operating normally.

Q: Can children or pets use the SunPlus EW300 Plus?

A: Yes, the device may be used for both children and pets when proper precautions are followed. Children should always use the device under adult supervision, and it should never be used on infants. Children over the age of 12 may generally follow the standard treatment guidelines for adults, while children under 12 should use approximately half the recommended treatment time.

Red light therapy is also commonly used for animals to help support issues such as discomfort, inflammation, and skin-related conditions. Before using the device on a pet, it is recommended to consult with a veterinarian to ensure treatment is safe and appropriate for your pet's specific needs.

Q: Should I use red light, near infrared light, or both together?

A: For most treatments, we recommend using both red light and near infrared (NIR) light together, as they work synergistically to support skin health, recovery, and

Welcome

Thank you for choosing All About Korean Beauty SunPlus EW300 Plus Red Light Therapy Device. This panel combines clinically recognized 660nm Red Light and 850nm Near Infrared Light to support recovery, skin wellness, circulation, relaxation, and overall wellness.

Designed for home and professional wellness environments, this device delivers targeted light therapy that may help support the body's natural recovery processes.

What's included

1. 1 x Protective Eyewear
2. 1 x US Power cord
3. 1 x User Manual
4. 1 x SunPlus EW300 Plus Red-Light Panel



Warnings

For safe and proper operation, read this manual carefully and follow all instructions. Store the device in a cool, dry location. To reduce the risk of injury or damage, DO NOT do the following:

- Do Not allow children to use the SunPlus EW300 Plus panel without adult supervision.
- Do Not use the device on individuals who are sleeping or unconscious.
- Do Not use the device together with products containing heat-producing ingredients, such as liniments, salves, ointments, creams, or balms, as this may cause burns or skin irritation.
- Do Not stare directly into the LEDs. Use protective eyewear if sensitive to light.
- Do Not Modify, disassemble, or attempt to repair the SunPlus EW300 Plus device. This device contains no parts that can be serviced by the user

To help ensure safe operation of the SunPlus EW300 Plus panel and reduce the risk of electrical hazards, please follow the safety guidelines below.

Do not perform any of the following actions while using or handling the SunPlus EW300 Plus:

- Do Not wash electrical components with water or other liquids.
- Do Not pull, carry, or lift the device using the power cord. If the power cord becomes damaged, discontinue use immediately.
- Do Not operate the device on or near heated surfaces.
- Do Not use the device in environments where flammable or combustible materials, gases, or vapors are present, as this may create a risk of fire or explosion.
- Do Not store the SunPlus EW300 Plus in temperatures exceeding 140°F (60°C).
- Do Not continue using the device if it appears damaged, as this could result in electrical shock or personal injury.
- Do Not use extension cords that are not rated to safely handle the device's power requirements.
- Do Not disconnect the device from the outlet by pulling on the cord. Always grip the plug itself when unplugging.
- Do Not block or cover the cooling fans while the device is in operation.
- Do Not place or store the device where it could fall into water or other liquids or expose it to running water. If the device comes into contact

make any medical claims or guarantees regarding the effectiveness of this product for physical, skin-related, or mental health conditions.

Warranty

All About Korean Beauty warrants that the SunPlus EW300 Plus will operate and perform according to its published specifications for a period of one (1) year from the date of delivery. During the first year of the warranty period, coverage includes parts and labor for eligible repairs, subject to certain exclusions. If warranty service is required, the customer is responsible for return shipping costs to All About Korean Beauty during the warranty period.













This warranty does not apply to devices that have been abused, misused, altered, improperly installed, damaged, or repaired with non-approved replacement parts. All About Korean Beauty is not responsible for any incidental, consequential, or indirect damages resulting from the use or failure of the product.

After the warranty period has expired, customers are responsible for all repair-related parts and labor costs. For any repair or replacement completed either during or after the warranty period, All About Korean Beauty reserves the right to provide replacement parts or units in colors that may differ from the original product.


To request warranty service or support, please contact:

info@allaboutkoreanbeauty.com

Recommended Usage Guidelines

RECOMMENDED USAGE GUIDELINES			
TREATMENT GOAL	RECOMMENDED DISTANCE	SESSION DURATION	RECOMMENDED FREQUENCY
 <p>SKIN CARE & REJUVENATION Support a healthy, radiant complexion and skin renewal.</p>	<p>18 inches</p> 	 <p>10-15 minutes</p>	 <p>3-5x weekly</p>
 <p>MUSCLE RECOVERY & WELLNESS Aid post-workout recovery, reduce soreness and support muscle performance.</p>	<p>6-12 inches</p> 	 <p>10-15 minutes</p>	 <p>3-5x weekly</p>
 <p>JOINT & DEEP TISSUE SUPPORT Promote joint comfort, improve mobility and support deep tissue health.</p>	<p>6-12 inches</p> 	 <p>15-20 minutes</p>	 <p>3-5x weekly</p>

TIP: Consistency is key for best results. Always begin with shorter sessions and adjust treatment time based on comfort and individual response.



with water or becomes submerged, immediately disconnect it from the power outlet. When using the device near water, it is recommended to use an outlet protected by a Ground Fault Circuit Interrupter (GFCI).

- Do Not touch the device, power plug, or switches with wet hands, or operate the device on wet or damp skin

Cautions

For safety, All About Korean Beauty recommends speaking with your healthcare provider before using this device if you are sensitive to light or are taking medications that may increase light sensitivity, including certain cold and allergy medications, pain relievers, or medications used to treat infections. Consultation with a healthcare provider is also recommended before using the SunPlus EW300 Plus panel if you are pregnant, have suspicious or cancerous skin lesions, have recently received steroid injections, or are using topical steroid treatments.

If you experience any discomfort, irritation, or concerns while using the device, discontinue use immediately and consult your healthcare provider. When the device is not in use, unplug it from the power outlet to help prevent tripping hazards and ensure safe operation.

Maintenance

The SunPlus EW300 Plus panels require minimal maintenance beyond routine care and regular cleaning. To clean the device, unplug it from the power source and gently wipe the surface with a damp cloth or paper towel. Avoid using harsh chemicals, abrasive cleaners, or bleach, as these may damage the panel. Always ensure the device is completely dry before reconnecting and operating it again.


For safety purposes, the SunPlus EW300 Plus panels are intended for indoor use only.

Disclaimer

Before using this product, especially if you have a medical condition or are under medical care, consult with a licensed healthcare professional to determine whether this product is appropriate for you. This product is not intended to diagnose, treat, cure, or prevent any disease or medical condition. All About Korean Beauty does not

Quick Start Guide

- Plug the device into a standard power outlet.
- Turn the device on using the control panel.
- Select brightness, pulse, and timer settings.
- Position the panel 6–18 inches from the treatment area.
- Select Treatment Time, Brightness, and Pulse Frequency. Use for 10–20 minutes per session.

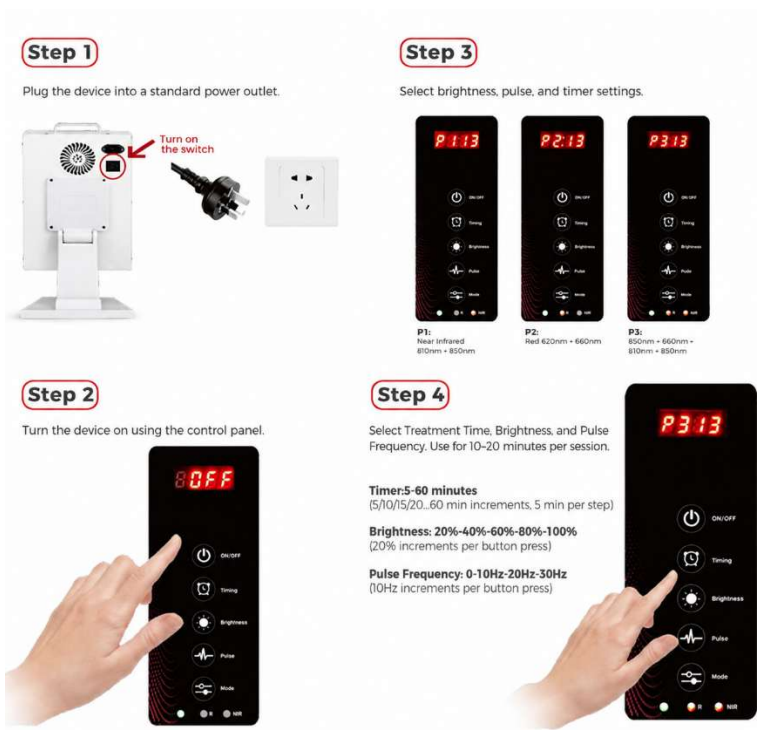

How Red-Light Therapy Works

Red light therapy (RLT) uses specific wavelengths of red and near-infrared light to stimulate cellular activity and promote healing. When the light penetrates the skin, it is absorbed by mitochondria, the “energy centers” of cells, boosting the production of adenosine triphosphate (ATP). This increased cellular energy helps repair tissues, reduce inflammation, and enhance circulation, while also promoting collagen and elastin production for healthier, more youthful-looking skin.

Key Benefits

- Muscle Recovery & Performance Support
- Skin Rejuvenation & Wellness
- Relaxation & Sleep Support
- Joint Comfort & Mobility
- Circulation & Recovery Support
- General Wellness & Energy Support

Technical Specifications

Light Panel Dimension:	12.3 x 9 x 2.6 inches (31 x 23 x 6.6cm)
Net Weight:	7.4lbs/3.34kg
Wavelength:	660nm, 850nm
LED Quantity:	60pcs * Dual Chip 5W
Function	Pulse Frequency: 0-10Hz-20Hz-30Hz Brightness: 20%-40%-60%-80%-100%